# **DIRECT SYSTEM**

VENEER OCCLUSIONVD COMPOSITE



beautiful innovation you can trust

# Bio-Esthetics and function in one appointment

edelweiss dentistry is a "think tank" driven by inspiration and technical know-how that conceptualizes and produces innovative systematic solutions that are developed together with partners within the dental industry and marketed via an exclusive dental distribution network. The edelweiss VENEER and OCCLUSIONVD SYSTEM are state of the art for modern and minimally invasive esthetic dentistry.

For the first time in the history of dentistry, it is now possible to work with prefabricated veneers made from Nano-Hybrid glass using modern laser technology. Never before has it been feasible to directly create the natural shape and youthful luminance of a tooth, so easily and perfectly in a single appointment. Its versatile application together with its time and cost saving procedure makes the edelweiss VENEER and OCCLUSIONVD SYSTEM a sound investment in the future with the best interest of the patient in mind. Convince yourself.



**Stephan Lampl**CEO, Founder & Inventor of edelweiss dentistry

Claudio Novelli
Clinical &
Scientific Director





# The Progressive Technical Concept

# edelweiss VENEER & OCCLUSIONVD

The laser-treated process combines the best of two worlds: a homogeneous, inorganic and maximum gloss surface with a sintered and dynamic Nano-Hybrid glass. Thus, an optimal integration in function and aesthetics is achieved.

This results in an anorganic biomechanical restoration in function and esthetics. The difference lies in its similarity with nature.

# **DIRECT SYSTEM**

VENEER OCCLUSIONVD COMPOSITE

#### **VENEER**

CERVICAL AREA: 0.2 mm

FACIAL AREA: 0.6 mm

INCISAL EDGE: 1.0 - 1.3 mm







### **OCCLUSIONVD**

MARGINAL AREA: 0.2 mm

OCCLUSION AREA: 0.6 - 1.3 mm



Solid and homogeneous glass phase

Laser sintered monobloc

	Technical Data	VENEER & OCCLUSIONVD
50 DEC 50	Flexural Strength	200 MPa
STATISTICS OF STATISTICS	Compressive Strength	550 MPa
	Flexural Modulus	20 GPa
Separate Sep	Surface Hardness	100 HV

(Source: internal data: R&D and QM department edelweiss dentistry)



# The Composition









# Step by Step VENEER



- 1 Non-prep preoperative situation with retraction cord
- **2** Etching with 37 % phosphoric acid
- **3** Applying adhesive to tooth surface



- 4 Adjusted & filled VENEER
- **5** Application of VENEER onto tooth
- 6 Position & press the VENEER into place, remove excess composite



- 7 Tack curing near the incisal edge with light curing unit
- 8 Cemented VENEER
- 9 Finishing of incisal edges


### **INDICATIONS**

- Anterior and posterior restorations
- Tooth discolorations
- Anatomical deformities
- Diastema
- Erosion & Attrition
- Customized Veneers, Crowns & Bridges
- Crowns & Crown facings
- Increasing the vertical dimension (CMD)

### **ADVANTAGES**

- One appointment
- Minimally invasive
- Natural look
- Long-lasting
- Biomechanical
- Biocompatible
- Antibacterial
- Economical

# Sizes & Shades

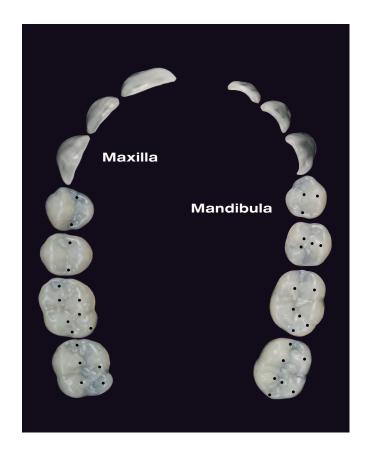



Claudio Novelli
P.H.D. D.M.D.
Clinical & Scientific Director

# The best for your patients

«By combining and adapting of the latest composite and laser technology, edelweiss dentistry has made a quantum leap into the future of minimally invasive dentistry. The laser process combines the best of two worlds: an homogeneous, inorganic and maximum gloss surface with a sintered and dynamic hybrid glass core. Thus, an optimal integration in function and esthetics is achieved. This, combined with maximum cost-effectiveness, has convinced me and my patients.»


«The difference lies in the high similarity with nature»

# Always the right size

The edelweiss VENEERs and OCCLU-SIONVDs are provided in different sizes. Using the edelweiss Sizing Guides, it is possible to easily determine the appropriate size of enamel shells.

VENEER Sizes: XS, S, M, L OCCLUSIONVD Sizes: S, M, L



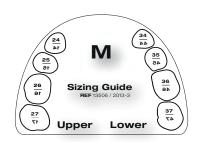

# **Shading**

edelweiss VENEERs and OCCLUSIONVD consist of the shade Enamel -

Vita Enamel AO. The respective dentin and enamel shades used to cement the edelweiss VENEER and OCCLUSIONVD will determine the final shade tone of the restoration.

Example: edelweiss VENEER (Vita Enamel Shade AO) cemented with Dentin Shade A3 (Vita Dentin Shade A3) will result in the Vita Dentin Shade A3.





ENAMEL
Shade of edelweiss
VENEER & OCCLUSIONVD



# Tailored to your workflow

# Direct or from the laboratory


The system from edelweiss dentistry offers the dentist three different treatment solutions which is tailored to the clinicians workflow.

Be your own creator and fit the VENEERs & OCCLUSIONVDs individually and directly on the patient. Tooth restorations in just one session can be realized simply and perfectly with the preshaped enamel shells. Or rely on our experience and get the VENEERS & OCCLUSIONVDs cutomized by our dental technicians.



# **Start Today**

Simply download our customization service lab sheet from the edelweiss website (www.edelweissdentistry.com/downloads). Fill out the form, attach it to your study model or impression and send it to our edelweiss laboratory. Perfect fitting was never that easy.

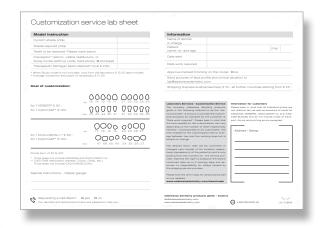



#### DIRECT

Adjust the VENEERs &
OCCLUSIONVDs directly on the
patient and then cement them using
adhesive technique as well as the
highly filled edelweiss NANO-HYBRID
COMPOSITE. The system allows
restorations within just one session.
This saves time and money for you
and your patient.



#### **SEMI-DIRECT**

You always have the possibility to adapt the VENEERs semi-directly by first grinding the VENEERs & OCCLUSIONVDs into shape, followed by bonding and then positioning with composite onto the unbonded tooth, removing surplus material and light curing. This allows you to prepare a perfect fit of all VENEERs before the final bonding and cementing.



# INDIRECT

You have no time to adjust the enamel shells directly or your case is quite complex? No problem.

Just send us the patients case related impression or study model together with the completed order form and get the VENEERs &

OCCLUSIONVDs or even crowns and bridges customized by our dental technicians.



# MALALIGNED TEETH / OLD RESTORATIONS

# edelweiss VENEER



**edelweiss VENEER** 13, 12, 11 & 21, 22, 23



# DISCOLORATION

# edelweiss VENEER



**edelweiss VENEER** 14, 13, 12, 11 & 21, 22, 23, 24



# DIASTEMA

# edelweiss VENEER



**edelweiss VENEER** 13, 12, 11 & 21, 22, 23



# ANATOMICAL DEFORMITIES / LENGTHENING OF THE INCISAL EDGE

# edelweiss VENEER



edelweiss VENEER
12 & 22
edelweiss N.H.-COMPOSITE
11 & 21, 13 & 23



# CROWN FACINGS / ABRASION

# edelweiss VENEER



**edelweiss VENEER** 14, 13, 12, 11 & 21, 22, 23, 24







- 1 Isolation with Opal Dam and etching with Ceramic Etch
- 2 Etching with 37 % phosphoric acid
- 3 Applying Bond



- 4 Tack curing near the incisal edge with light curing unit
- 5 Covering metal area with edelweiss EFFECT SHADE Opaque White
- 6 Application of VENEER
- 7 Position & press the VENEER into place

# FULL MOUTH REHABILITATION

# edelweiss VENEER & OCCLUSIONVD





# SINGLE RESTAURATIONS WITH OCCLUSIONVD

# edelweiss OCCLUSIONVD


# Direct ceramic crown facing



# Direct replacement of large amalgam



# Direct post endodontic tooth onlay

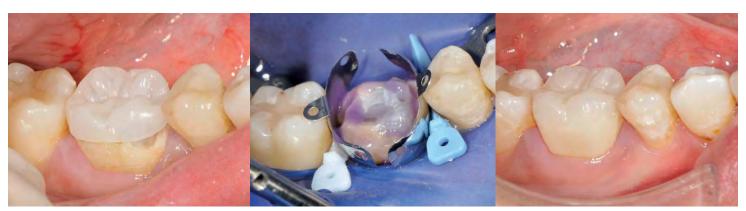

# Expertise & Knowledge

# Live-Workshops for your success

Experience the edelweiss system live on one of our Workshops with lectures and hands-on demonstrations! Get in touch with your sales representative or go to

www.edelweissdentistry.com/ education-center





Get your edelweiss compendium and other advertising materials as a PDF download and show your patients the versatile options provided by the edelweiss system.



# Online education on our edelweiss website



Visit our edelweiss video library!

Explore our step by step videos and quick tips on the edelweiss website and learn how to optimally treat your patients within a short time.

www.edelweissdentistry.com/videos



Marc Bachmann Webinar DT Study Club Bio-Esthetic & Function in one appointment



Sanette Muller & Stephan Lampl Amelogenesis Imperfecta



Stephan Lampl & Claudio Novelli Full Mouth Rehabilitation Full Mouth Rehabilitation Full Mouth Rehabilitation Characterization of with edelweiss VENEERs My story / & OCCLUSIONVDs



Stephan Lampl & Claudio Novelli Robert Lucero



**Didier Dietschi** Implementation & edelweiss VENEERs

# Step by Step OCCLUSIONVE



Preoperative situation. Occlusal view.



Prepared and roughened crown.



Roughen inner surface of the OCCLUSIONVD.



Insertion of the OCCLUSIONVD.



Positioning and fit check of the OCCLUSIONVD in final occlusal position.



Positioning and fit check of established position in occlusion.



Removal of the semi direct fitted OCCLUSIONVD.



Basic view of the semi direct fitted OCCLUSIONVD.



Isolation with Wedge Guards and etching for 30 sec.



Insertion and final positioning of OCCLUSIONVD.



Final fit check of OCCLUSIONVD.



Adapt margins. Removal of surplus material.



Finishing interproximal areas using polishing strip.




Visual observation.



Reocclusion.



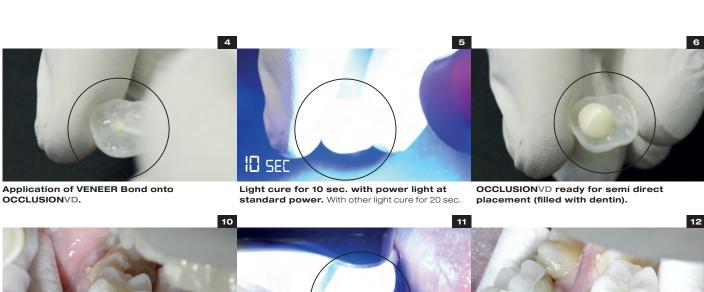



**Light cure with power light for 10 sec.** With other light cure 20 sec.





Placing a puddle coat of Bond on crown.



With power light cure for 10 sec. at standard Adapt Nano Hybrid Composite into power. With other light cure for 20 sec.



OCCLUSIONVD.



With power light cure for 10 sec. at standard Removal of Wedge Guards. power. With other light cure for 20 sec.





Finish cervical and interproximal margins using Soflex Discs.



Reocclusion.



For high shine polish OCCLUSIONVD margins with KENDA polishers.



Final postoperative situation.



# ADVANCED DIRECT SYSTEM VENEER & OCCLUSIONVD TOOLBOX

REF 13700	
36	edelweiss VENEERs
1	VENEER SET S up (1 x 14, 13, 12, 11, 21, 22, 23, 24 S)
1	VENEER SET M up (1 x 14, 13, 12, 11, 21, 22, 23, 24 M)
1	VENEER SET L up (1 x 14, 13, 12, 11, 21, 22, 23, 24 L)
1	VENEER SET S low (1 × 43, 42, 41, 31, 32, 33 S)
1	VENEER SET M low (1 x 43, 42, 41, 31, 32, 33 M)
48	edelweiss OCCLUSIONVD
2	OCCLUSIONVD SETS up (1 x 17, 16, 15, 14 S) (1 x 24, 25, 26, 27 S)
2	OCCLUSIONVD SET M up $(1 \times 17, 16, 15, 14 \text{ M})$ $(1 \times 24, 25, 26, 27 \text{ M})$
2	OCCLUSIONVD SET L up (1 x 17, 16, 15, 14 L) (1 x 24, 25, 26, 27 L)
2	OCCLUSIONVD SETS low (1 x 47, 46, 45, 44 S) (1 x 34, 35, 36, 37 S)
2	OCCLUSIONVD SET M low (1 $\times$ 47, 46, 45, 44 M) (1 $\times$ 34, 35, 36, 37 M)
2	OCCLUSIONVD SET L low (1 x 47, 46, 45, 44 L) (1 x 34, 35, 36, 37 L)
25	edelweiss NANO-HYBRID COMPOSITE DENTIN SHADE
5	Dentin A0 0.3 g Tip
5	Dentin A1 0.3 g Tip
5	Dentin A2 0.3 g Tip
5	Dentin A3 0.3 g Tip
5	Dentin A3.5 0.3 g Tip
5	edelweiss NANO-HYBRID COMPOSITE ENAMEL SHADE
5	Enamel 0.3 g Tip
1	edelweiss NANO-HYBRID COMPOSITE FLOW
1	Enamel Flowable 1.5 g Syringe
2	edelweiss EFFECT SHADES

- Opaque White 1.5 g Syringe
- Effect Blue 1.5 g Syringe

# edelweiss BOND

VENEER Bond 5 ml Bottle

### ACCESSORIES

Dentin & Enamel Shade Guide VENEER & OCCLUSIONVD Sizing Guide up & low FLOW Application Tips (10x) Step by Step Booklet Instructions for use



Removable cover for more flexibility





# STARTER DIRECT SYSTEM VENEER TOOLBOX REF 13710

20	edelweiss VENEERs
1	VENEER SET S up (1 x 13, 12, 11, 21, 22, 23 S)
1	VENEER SET M up (1 x 14, 13, 12, 11, 21, 22, 23, 24 M)
1	VENEER SET L up (1 x 13, 12, 11, 21, 22, 23 L)
25	edelweiss NANO-HYBRID COMPOSITE DENTIN SHADE

5	Dentin A0 0.3 g Tip
5	Dentin A1 0.3 g Tip
5	Dentin A2 0.3 g Tip
5	Dentin A3 0.3 g Tip

Dentin A3.5 0.3 g Tip

# 5 edelweiss NANO-HYBRID COMPOSITE ENAMEL SHADE

5 Enamel 0.3 g Tip

5

# 1 edelweiss NANO-HYBRID COMPOSITE FLOW

1 Enamel Flowable 1.5 g Syringe

### 2 edelweiss EFFECT SHADES

Opaque White 1.5 g Syringe
 Effect Blue 1.5 g Syringe

### 1 edelweiss BOND

1 VENEER Bond 5 ml Bottle

#### **ACCESSORIES**

Dentin & Enamel Shade Guide
VENEER Sizing Guide up & low
FLOW Application Tips (10x)
Step by Step Booklet
Instructions for use



# STARTER DIRECT SYSTEM OCCLUSIONVD TOOLBOX REF 13720

32	edelweiss OCCLUSIONVD
2	OCCLUSIONVD SET S up (1 x 17, 16, 15, 14 S) (1 x 24, 25, 26, 27 S)
2	OCCLUSION√D SET M up (1 x 17, 16, 15, 14 M) (1 x 24, 25, 26, 27 M)
2	OCCLUSIONVD SET S low $(1 \times 47, 46, 45, 44  S) (1 \times 34, 35, 36, 37  S)$
2	OCCLUSIONVD SET M low (1 x 47, 46, 45, 44 M) (1 x 34, 35, 36, 37 M)
25	edelweiss NANO-HYBRID COMPOSITE DENTIN SHADE
5	Dentin A0 0.3 g Tip
_	

	DENTIN SHADE
5	Dentin A0 0.3 g Tip
5	Dentin A1 0.3 g Tip
5	Dentin A2 0.3 g Tip
5	Dentin A3 0.3 g Tip
5	Dentin A3.5 0.3 g Tip

# 5 edelweiss NANO-HYBRID COMPOSITE ENAMEL SHADE

5 Enamel 0.3 g Tip

# 1 edelweiss NANO-HYBRID COMPOSITE FLOW

1 Enamel Flowable 1.5 g Syringe

#### 1 edelweiss BOND

1 VENEER Bond 5 ml Bottle

#### **ACCESSORIES**

Dentin & Enamel Shade Guide

OCCLUSIONVD Sizing Guide up & Iow

FLOW Application Tips (5x)

Step by Step Booklet

Instructions for use

# Refills

	ENEER SET REFILLS	REF 10161	edelweiss VENEER 21 L
REF 13000	edelweiss VENEER SET 14-24 S up (1 x 14, 13, 12, 11, 21, 22, 23, 24 S)	REF 14128	edelweiss VENEER 22 XS
DEE 10010	,	REF 10162	edelweiss VENEER 22 S
REF 13010	edelweiss VENEER SET 14-24 M up (1 x 14, 13, 12, 11, 21, 22, 23, 24 M)	REF 10163	edelweiss VENEER 22 M
REF 13020	edelweiss VENEER SET 14-24 L up	REF 10164	edelweiss VENEER 22 L
TIET 10020	(1 x 14, 13, 12, 11, 21, 22, 23, 24 L)	REF 14129	edelweiss VENEER 23 XS
REF 13030	edelweiss VENEER SET 44-34 S low	REF 10165	edelweiss VENEER 23 S
	(1 × 44, 43, 42, 41, 31, 32, 33, 34 S)	REF 10166	edelweiss VENEER 23 M
REF 13040	edelweiss VENEER SET 44-34 M low	REF 10167	edelweiss VENEER 23 L
	(1 × 44, 43, 42, 41, 31, 32, 33, 34 M)	REF 13106	edelweiss VENEER 24 S
REF 14121	edelweiss VENEER SET 13-23 XS up	REF 13107	edelweiss VENEER 24 M
	(1 x 13, 12, 11, 21, 22, 23 XS)	REF 13108	edelweiss VENEER 24 L
REF 10100	edelweiss VENEER SET 13-23 S up (1 x 13, 12, 11, 21, 22, 23 S)	REF 13109	edelweiss VENEER 25 S
DEE 10110		REF 13110	edelweiss VENEER 25 M
REF 10110	edelweiss VENEER SET 13-23 M up (1 x 13, 12, 11, 21, 22, 23, M)	REF 13111	edelweiss VENEER 25 L
REF 10120	edelweiss VENEER SET 13-23 L up	REF 10168	edelweiss VENEER 31 S
1121 10120	(1 × 13, 12, 11, 21, 22, 23 L)	REF 10169	edelweiss VENEER 31 M
REF 10130	edelweiss VENEER SET 43-33 S low	REF 10170	edelweiss VENEER 32 S
	(1 x 43, 42, 41, 31, 32, 33 S)	REF 10171	edelweiss VENEER 32 M
REF 10140	edelweiss VENEER SET 43-33 M low	REF 10172	edelweiss VENEER 33 S
	(1 x 43, 42, 41, 31, 32, 33 M)	REF 10173	edelweiss VENEER 33 M
		REF 13112	edelweiss VENEER 34 S
odobycios V	ENEER SINGLE REFILLS	REF 13113	edelweiss VENEER 34 M
		REF 13114	edelweiss VENEER 35 S
REF 14122	edelweiss VENEER 11/21 XS (1 x 11, 21 XS)	REF 13115	edelweiss VENEER 35 M
REF 10180	edelweiss VENEER 11/21 S (1 x 11, 21 S)	REF 10174	edelweiss VENEER 41 S
REF 10181	edelweiss VENEER 11/21 M (1 x 11, 21 M)	REF 10175	edelweiss VENEER 41 M
REF 10182	edelweiss VENEER 11/21 L (1 x 11, 21 L)	REF 10176	edelweiss VENEER 42 S
REF 14123	edelweiss VENEER 12/22 XS (1 x 12, 22 XS)	REF 10177	edelweiss VENEER 42 M
REF 10183	edelweiss VENEER 12/22 S (1 x 12, 22 S)	REF 10178	edelweiss VENEER 43 S
REF 10184	edelweiss VENEER 12/22 M (1 x 12, 22 M)	REF 10179	edelweiss VENEER 43 M
REF 10185	edelweiss VENEER 12/22 L (1 x 12, 22 L) edelweiss VENEER 11 XS	REF 13116	edelweiss VENEER 44 S
REF 14124 REF 10150	edelweiss VENEER 11 S	REF 13117	edelweiss VENEER 44 M
REF 10150		REF 13118	edelweiss VENEER 45 S
REF 10151	edelweiss VENEER 11 M edelweiss VENEER 11 L	REF 13119	edelweiss VENEER 45 M
	edelweiss VENEER 11 L edelweiss VENEER 12 XS		
REF 14125			
REF 10153	edelweiss VENEER 12 S	edelweiss (	OCCLUSIONVD SET REFILLS
REF 10154 REF 10155	edelweiss VENEER 12 M	REF 13170	OCCLUSIONVD SET 17-27 S up
	edelweiss VENEER 12 L		(1 × 17, 16, 15, 14, 24, 25, 26, 27 S)
REF 14126	edelweiss VENEER 13 XS edelweiss VENEER 13 S	REF 13175	OCCLUSIONVD SET 17-27 M up
REF 10156			(1 x 17, 16, 15, 14, 24, 25, 26, 27 M)
REF 10157	edelweiss VENEER 13 M	REF 13180	OCCLUSIONVD SET 17-27 L up
REF 10158	edelweiss VENEER 13 L	DEE 10105	(1 x 17, 16, 15, 14, 24, 25, 26, 27 L)
REF 13100	edelweiss VENEER 14 S	REF 13185	OCCLUSIONVD SET 47-37 S low (1 x 47, 46, 45, 44, 34, 35, 36, 37 S)
REF 13101	edelweiss VENEER 14 M	REF 13190	OCCLUSIONVD SET 47-37 M low
REF 13102	edelweiss VENEER 14 L	1161 10190	(1 x 47, 46, 45, 44, 34, 35, 36, 37 M)
REF 13103	edelweiss VENEER 15 S	REF 13195	OCCLUSIONVD SET 47-37 L low
REF 13104	edelweiss VENEER 15 M		(1 × 47, 46, 45, 44, 34, 35, 36, 37 L)
REF 13105	edelweiss VENEER 15 L	REF 13150	OCCLUSIONVD SET 17-14 S up
REF 14127	edelweiss VENEER 21 XS		(1 × 17, 16, 15, 14 S)
REF 10159	edelweiss VENEER 21 S	REF 13151	OCCLUSIONVD SET 17-14 M up
REF 10160	edelweiss VENEER 21 M		(1 × 17, 16, 15, 14 M)



# **DIRECT SYSTEM**

REF 13152	OCCLUSIONVD SET	17-14 L up	REF 13232	OCCLUSIONVD	36 L
DEE 10150	(1 × 17, 16, 15, 14 L)	0.4.07.0	REF 13233	OCCLUSIONVD	37 S
REF 13153	OCCLUSIONVD SET (1 x 24, 25, 26, 27 S)	24-27 S up	REF 13234	OCCLUSIONVD	37 M
REF 13154	OCCLUSIONVD SET	24-27 M up	REF 13235	OCCLUSIONVD	37 L
	(1 x 24, 25, 26, 27 M)		REF 13236	OCCLUSIONVD	44 S
REF 13155	EF 13155 OCCLUSIONVD SET 2	1-27 L up	REF 13237	OCCLUSIONVD	44 M
	(1 × 24, 25, 26, 27 L)		REF 13238	OCCLUSIONVD	44 L
REF 13156	OCCLUSIONVD SET	34-37 S low	REF 13239	OCCLUSIONVD	45 S
	(1 × 34, 35, 36, 37 S)		REF 13240	OCCLUSIONVD	45 M
REF 13157	OCCLUSIONVD SET (1 x 34, 35, 36, 37 M)		REF 13241	OCCLUSIONVD	45 L
DEE 10150			REF 13242	OCCLUSIONVD	46 S
REF 13158	OCCLUSIONVD SET (1 x 34, 35, 36, 37 L)	34-37 L IOW	REF 13243	OCCLUSIONVD	46 M
REF 13159		47-44 S low	REF 13244	OCCLUSIONVD	46 L
1121 10100	(1 × 47, 46, 45, 44 S)	77 1101000	REF 13245	OCCLUSIONVD	47 S
REF 13160	OCCLUSIONVD SET	47-44 M low	REF 13246	OCCLUSIONVD	47 M
	(1 × 47, 46, 45, 44 M)		REF 13247	OCCLUSIONVD	47 L
REF 13161	OCCLUSIONVD SET (1 x 47, 46, 45, 44 L)	47-44 L low			
	CCLUSIONVD SING	_	_	RID COMPOSITE DE	
REF 13200	OCCLUSIONVD	14 S	REF 10201	Dentin A0	15 x 0.3 g Tip
REF 13201	OCCLUSIONVD	14 M	REF 10211	Dentin A1	15 x 0.3 g Tip
REF 13202	OCCLUSIONVD	14 L	REF 10221	Dentin A2	15 x 0.3 g Tip
REF 13203	OCCLUSIONVD	15 S	REF 10231	Dentin A3	15 x 0.3 g Tip
REF 13204	OCCLUSIONVD	15 M	REF 10241	Dentin A3.5	15 x 0.3 g Tip
REF 13205	OCCLUSIONVD	15 L			
REF 13206	OCCLUSIONVD	16 S	NANO-HYBF	RID COMPOSITE EN	AMEL SHADE
REF 13207	OCCLUSIONVD	16 M	REF 10301	Enamel	10 x 0.3 g Tip
REF 13208	OCCLUSIONVD	16 L			3 1
REF 13209	OCCLUSIONVD	17 S			
REF 13210	OCCLUSIONVD	17 M	EFFECT SHA	_	45 0 1
REF 13211	OCCLUSIONVD	17 L	REF 13760	Effect Blue	1.5 g Syringe
REF 13212	OCCLUSIONVD	24 S	REF 13750	Opaque White	1.5 g Syringe
REF 13213	OCCLUSIONVD	24 M	REF 14702	Effect Ice	1.5 g Syringe
REF 13214	OCCLUSIONVD	24 L			
REF 13215	OCCLUSIONVD	25 S	NANO-HYBE	RID COMPOSITE FLO	ow
REF 13216	OCCLUSIONVD	25 M	REF 13770	Enamel Flowable	1.5 g Syringe
REF 13217	OCCLUSIONVD	25 L			
REF 13218	OCCLUSIONVD	26 S	BONDING		
REF 13219	OCCLUSIONVD	26 M	REF 10521	VENEER Bond	5 ml Bottle
REF 13220	OCCLUSIONVD	26 L	110021	V EI VEEL I DONG	O IIII DOLLIO
REF 13221	OCCLUSIONVD	27 S			
REF 13222	OCCLUSIONVD	27 M	ACCESSORI		
REF 13223	OCCLUSIONVD	27 L	REF 14600	Dentin & Enamel Sha	
REF 13224	OCCLUSIONVD	34 S	REF 13508	VENEER Sizing Guide up & low	

REF 13225 OCCLUSIONVD

REF 13226 OCCLUSIONVD

REF 13228 OCCLUSIONVD

REF 13229 OCCLUSIONVD

REF 13231 OCCLUSIONVD

REF 13230 OCCLUSIONVD

REF 13227 OCCLUSIONVD

34 M

34 L

35 S

35 M

35 L

36 S

36 M

# 19

REF 13509 OCCLUSIONVD Sizing Guide up & low

REF 814027 Step by Step Booklet

FLOW Application Tips (10x)

REF 10621

PROMOTION TOOL

Available upon request!



beautiful innovation you can trust

# 816029/ 2022-05



edelweiss dentistry ® and the edelweiss Logo ® are registered trademarks of edelweiss dentistry products gmbh • Austria

© 2022 edelweiss dentistry products gmbh • Austria. All material contained in this brochure are the property of edelweiss dentistry products gmbh • Austria. No part of this publication may be reproduced, distributed, or transmitted in any form or by any means, without the prior written permission of the edelweiss dentistry products gmbh • Austria. All Rights Reserved.